



Yes, you can.®

➤ Invacare® **Birdie®** and **Birdie® Compact**

Transportability and patient comfort go hand-in-hand

The innovative boom design of the *Invacare® Birdie®* and the *Birdie® Compact* offers a truly comfortable lift and transfer to or from bed, chair or even the floor. The lifters are designed for maximum convenience and offer a folding format that can be carried out easily and without the need for tools. In addition, the lifter can be dismantled into two parts if needed – also without tools.

3 year warranty





➤ **Total user comfort**

Allowing a 360° range of movement even when the hoist is in the highest position, the **Birdie**® provides a very comfortable lift and transfer. This allows the patient to be easily rotated in the hoist which allows placement into the bed or wheelchair much easier for the carer. Optimising space for the client creates a high level of comfort and ensures that the clients' knees are not too close to the actuator of the hoist, which means that there is a very low risk of injury.

➤ **Ease of Handling: simple to fold, unfold and dismantle**

The hoist can easily be folded and unfolded without the need for tools. The hoist will require minimal space when stored and is very easy to push and transport*. If required the mast and boom can be easily detached from the base (no tools required).

➤ **TÜV certification for maximum security**

For total security, the hoists have a snap lock system to minimise the risk of accidents. The hoists are CE certified in line with the Medical Device Directive 93/42/ECC. To ensure the hoists comply with EN ISO 10535, the hoists have been rigorously tested by Invacare and the German TÜV institute, an accredited external body.



➤ **Configured to client needs**

Clients can enjoy maximum flexibility with a range of spreader bars and a very comfortable sling range which allows for individual adjustment of the transfer.

* For transportation between external locations, Invacare® recommend the **Birdie**® and **Birdie**® **Compact** be dismantled into two parts.



Features



Long boom – high user comfort

The client can easily be rotated and positioned. Reduces risk of injuries to the lower limbs.



Easy to fold and dismantle

Ease of transport and storage. No tools required when folding / unfolding or when dismantling the hoist.



Snap lock

Lock snaps into place. A magnetic pin ensures optimal security.



Emergency stop, emergency lowering and integrated charger

Manual or electric emergency lowering (compact version features only electric lowering). To charge the hoist, simply plug in.



Manual leg spread

Simply step on the pedal to adjust the leg spread (**Birdie**® and **Birdie Compact**).



Clip on spreader bar

Allows the spreader bar to be changed quickly.



Options and Accessories



2 point quick release spreader bars
Easy to attach and detach. Available widths: 45 cm, 55 cm (35 mm also available).



4 point quick release spreader bars
Easy to attach and detach. Available widths: 45 cm and 55 cm.



Electrical leg spread
Available for *Invacare® Birdie®*.



Scale
Approved by the medical standard EN 45501. Certified and sealed.



Handle for leg spread
Handle can easily be attached to operate the manual leg spread.





Product Part Number

3000500.E2A00-0125
3000500.E2C00-0125
3000500.D2C00-0125
3000500.D2A00-0125
3000510.D2C00-0125
3000510.D2A00-0125


Product Description

Birdie® with electrical leg opening and 2 point spreader bar
Birdie® with electrical leg opening and 4 point spreader bar
Birdie® with manual leg opening and 4 point spreader bar
Birdie® with manual leg opening and 2 point spreader bar
Birdie® Compact with manual leg opening and 4 point spreader bar
Birdie® Compact with manual leg opening and 2 point spreader bar

Technical data

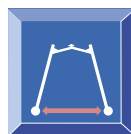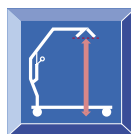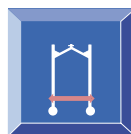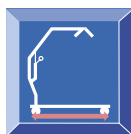

| | | | | | | |
|------------------------|--------|-------|-------|-----------------|---------|--------|
| Birdie® | 115 cm | 68 cm | 42 kg | 46.5 - 173.5 cm | 11.5 cm | 104 cm |
| Birdie® Compact | 100 cm | 55 cm | 36 kg | 54 - 161 cm | 11.5 cm | 89 cm |

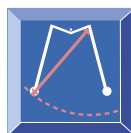

| | | | | |
|------------------------|--------|-----------|--------|--------|
| Birdie® | 170 kg | 40 cycles | 2,9 Ah | 140 cm |
| Birdie® Compact | 150 kg | 40 cycles | 2,9 Ah | 107 cm |

Colour of frame



Light grey / Dark grey

Birdie · UK · 04/09 · Art. No. 1504431

© Invacare® International Srl 2009. All rights reserved.

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Limited
South Road · Bridgend Industrial Estate · BRIDGEND CF31 3PY
Tel. +44 (0) 1656 647327 · Fax +44 (0) 1656 649016
E-mail: uk@invacare.com · www.invacare.co.uk
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 · Seatown Business Campus
Seatown Road · SWORDS · County Dublin
Tel. (353) 181 070 84 · Fax (353) 181 070 85
Email: Ireland@invacare.com · www.invacare.ie